

Cover Screws

Instructions for use



Important: Please read.

Disclaimer of liability:

This product is part of an overall concept and may only be used in conjunction with the associated original products according to the instructions and recommendation of Nobel Biocare. Non-recommended use of products made by third parties in conjunction with Nobel Biocare products will void any warranty or other obligation, express or implied, of Nobel Biocare. The user of Nobel Biocare products has the duty to determine whether or not any product is suitable for the particular patient and circumstances. Nobel Biocare disclaims any liability, express or implied, and shall have no responsibility for any direct, indirect, punitive or other damages, arising out of or in connection with any errors in professional judgment or practice in the use of Nobel Biocare products. The user is also obliged to study the latest developments in regard to this Nobel Biocare product and its applications regularly. In cases of doubt, the user has to contact Nobel Biocare. Since the utilization of this product is under the control of the user, they are his/her responsibility. Nobel Biocare does not assume any liability whatsoever for damage arising thereof. Please note that some products detailed in this Instruction for Use may not be regulatory cleared, released or licensed for sale in all markets.

Description:

A premanufactured device intended to be directly connected to the endosseous dental implant.

Internal conical connection for: NobelActive®, NobelReplace® CC and NobelParallel™ CC. Internal tri-channel connection for: NobelReplace®, Replace Select™, NobelSpeedy® Replace, NobelReplace® Platform Shift.

External hex connection for: Brånemark System®, NobelZygoma™, Zygoma Implant and NobelSpeedy® Groovy.

Intended use:

Dental implant cover screws are to be used in the upper or lower jaw connected to the endosseous implant to protect the internal threads and implant head during the healing phase.

Indication:

To be used when applicable together with the implant during healing in order to protect the implant platform and internal threads from overgrowth of bone.

Contraindications:

Cover Screw is contraindicated for patients:

- who are medically unfit for an oral surgical procedure.
- in whom adequate sizes, numbers or desirable positions of implants are not reachable to achieve safe support of functional or eventually parafunctional loads.
- allergic or hypersensitive to commercially pure titanium or titanium alloy Ti-6Al-4V (titanium, aluminum, vanadium).

Cautions:

Pre-operative hard tissue or soft tissue deficits may yield a compromised esthetic result or unfavorable implant angulations.

To secure the long term treatment outcome it is advised to provide comprehensive regular patient follow up after implant treatment and to inform about appropriate oral hygiene. All instruments and tooling used in surgery must be maintained in good condition and care must be taken that instrumentation does not damage implants or other components. Special attention has to be given to patients who have local or systemic factors that could interfere with the healing process of either bone or soft tissue or the osseointegration process (e.g., cigarette smoking, poor oral hygiene, uncontrolled diabetes, oro-facial radiotherapy, steroid therapy, infections in the neighboring bone).

Special caution is advised in patients who receive bisphosphonate therapy.

In general, implant placement and prosthetic design must accommodate individual patient conditions. In case of bruxism or unfavorable jaw relationships reappraisal of the treatment option may be considered.

Because of the small size of the devices, care must be taken that they are not swallowed or aspirated by the patient.

It is strongly recommended that clinicians, new as well as experienced implant users, always go through special training before undertaking a new treatment method. Nobel Biocare offers a wide range of courses for various levels of knowledge and experience. For more info please visit www.nobelbiocare.com.

Working the first time with a colleague, experienced with the new device/treatment method, avoids eventual complications. Nobel Biocare has a global network of mentors available for this purpose.


Close cooperation between surgeon, restorative dentist and dental laboratory technician is essential for a successful implant treatment.

Handling instructions:

Clinical procedure:

1. Select appropriate cover screw.
2. Connect to implant and hand tighten using Unigrip™ Screwdriver (A).

Note: For Brånemark System® cover screw a hex screwdriver should be used for tightening.

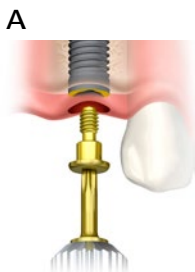

Materials:

Cover Screw for implants with External Hex connection: Commercially pure titanium.

Cover Screw for implants with Internal conical connection, Internal tri-channel connection and Zygoma implants: Titanium alloy 90% Ti, 6% Al, 4%V.

Cleaning and sterilization instructions:

Cover Screw is delivered sterile and for single use only prior to the labeled expiration date.

Warning: Do not use device if the packaging has been damaged or previously opened.

Caution: Cover Screw is a single use product not intended to be reprocessed. Reprocessing could cause loss of mechanical, chemical and/or biological characteristics. Reuse could cause cross contamination.

MR safety information:

Note: Only the Conical Connection Wide Platform Cover Screw has been assessed as MR Conditional. The other platforms and sizes have not been evaluated for safety and compatibility in the MR environment and have not been tested for heating or migration in the MR environment.

Non-clinical testing has demonstrated that the product is MR conditional. A patient with this device can be safely scanned in a MR system meeting the following conditions:

- Static magnetic field of 1.5-Tesla and 3-Tesla only.
- Maximum spatial gradient magnetic field of 4000 Gauss/cm or less (40 T/m).
- Maximum MR system reported, whole body averaged specific absorption rate (SAR) of 2 W/kg (Normal Operating Mode) or 4 W/kg (First Level Controlled Mode).

Under the scan conditions defined above, the product is expected to produce a maximum temperature rise of 4.1° C after 15 minutes of continuous scanning.

In non-clinical testing, the image artifact caused by the device extends approximately 30mm from the product when imaged with a gradient echo pulse sequence a 3 Tesla MRI system.

Should there be no MR symbol on the product label, please note that the product has not been evaluated for safety and compatibility in the MR environment. The product has not been tested for heating or migration in the MR environment.

For additional information on Cleaning and Sterilization and Magnetic Resonance Imaging, please consult the "Cleaning & Sterilization Guidelines for Nobel Biocare Products including MRI Information" available at www.nobelbiocare.com/sterilization or request latest printed version from a Nobel Biocare representative.

Storage and handling:

The product must be stored in a dry place in the original packaging at room temperature and not exposed to direct sunlight. Incorrect storage may influence device characteristics leading to failure.

Disposal:

Disposal of the device shall follow local regulations and environmental requirements, taking different contamination levels into account.




Manufacturer: Nobel Biocare AB, Box 5190, 402 26
Västra Hamngatan 1, 411 17 Göteborg, Sweden.

Phone: +46 31 81 88 00. Fax: +46 31 16 31 52. www.nobelbiocare.com

Canada license exemption: Please note that not all products may have been licensed in accordance with Canadian law.

Prescription device: Rx only

Caution: Federal law restricts this device to sale by or on the order of a licensed physician or dentist.

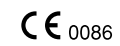

Rx Only



Magnetic
resonance
conditional



Sterile using
irradiation



Consult
instructions
for use



Use-by
date



Do not
re-use



Batch code



Do not use
if package is
damaged

EN All rights reserved.

Nobel Biocare, the Nobel Biocare logotype and all other trademarks used in this document are, if nothing else is stated or is evident from the context in a certain case, trademarks of Nobel Biocare. Product images in this folder are not necessarily to scale.