

Cover Screws



Important – Disclaimer of Liability

This product is part of an overall concept and may only be used in conjunction with the associated original products according to the instructions and recommendation of Nobel Biocare. Non-recommended use of products made by third parties in conjunction with Nobel Biocare products will void any warranty or other obligation, express or implied, of Nobel Biocare. The user of Nobel Biocare products has the duty to determine whether or not any product is suitable for the particular patient and circumstances. Nobel Biocare disclaims any liability, express or implied, and shall have no responsibility for any direct, indirect, punitive or other damages, arising out of or in connection with any errors in professional judgment or practice in the use of Nobel Biocare products. The user is also obliged to study the latest developments in regard to this Nobel Biocare product and its applications regularly. In cases of doubt, the user has to contact Nobel Biocare. Since the utilization of this product is under the control of the user, they are his/her responsibility. Nobel Biocare does not assume any liability whatsoever for damage arising thereof.

Please note that some products detailed in this Instruction for Use may not be regulatory cleared, released or licensed for sale in all markets.

Description

Cover screws are a component that covers the implant platform and prevents tissue overgrowth during the healing phase of the implant. The threaded portion of the Cover Screw fits inside the internal thread of the implant, while the head of the Cover Screw covers the top surface of the implant (the implant head).

The Cover Screw Nobel Biocare N1™ TCC features two parts: the main body and the internal screw.

Nobel Biocare products are intended and available to be used in a variety of configurations. For further information refer to Nobel Biocare publication Compatibility Information by navigating to ifu.nobelbiocare.com.

Intended Use/Intended Purpose

Dental Implant Cover Screws

Dental implant Cover Screws are to be used in the upper or lower jaw connected to the endosseous implant to protect the internal threads and implant head during the healing phase.

Cover screw Nobel Biocare N1™

Intended to be temporarily connected to an endosseous dental implant to protect the implant connection interface during bone healing.

Indications

Dental Implant Cover Screws

To be used when applicable together with the implant during healing in order to protect the implant platform and internal threads from overgrowth of bone.

Cover Screws Nobel Biocare N1™

Cover screws are indicated for use with implants placed in the maxilla, mandible, or zygomatic bone, as per the indications for the respective implant system.

Contraindications

It is contraindicated to use Cover Screws in:

- Patients who are medically unfit for an oral surgical procedure.
- Patients in whom adequate sizes, numbers, or desirable positions of implants are not reachable to achieve safe support of static and dynamic loads.
- Patients allergic or hypersensitive to commercially pure titanium or titanium alloy Ti-6Al-4V (titanium, aluminum, vanadium) and DLC (Diamond Like Carbon) coating.

For contraindications specific to the implants, refer to the Nobel Biocare Instructions for Use for the Implant Systems.

Materials

Cover Screw for implants with Internal Conical Connection, Internal Tri-Channel and External Hex, connection

Screw: Titanium alloy composed of Titanium-6 Aluminium-4 Vanadium ELI (Extra Low Interstitial) Alloy. Detailed chemical composition of the Titanium alloy is Titanium balanced with 5.50 - 6.50 wt.% Aluminium, 3.5-4.5 wt.% Vanadium, max. 0.25 wt.% Iron, max. 0.13 wt.% Oxygen, max. 0.08 wt.% Carbon, max. 0.05 wt.% Nitrogen, max. 0.012 wt.% Hydrogen, (max. -maximum value).

Cover Screws for the Nobel Biocare N1 Implant System

Screw: Titanium alloy composed of Titanium-6 Aluminium-4 Vanadium ELI (Extra Low Interstitial) Alloy. Detailed chemical composition of the Titanium alloy is Titanium balanced with 5.50 - 6.50 wt.% Aluminium, 3.5-4.5 wt.% Vanadium, max. 0.25 wt.% Iron, max. 0.13 wt.% Oxygen, max. 0.08 wt.% Carbon, max. 0.05 wt.% Nitrogen, max. 0.012 wt.% Hydrogen, (max. -maximum value). Screw is partly coated by Diamond like Carbon Coating. The coating is metal containing carbon coating, containing tungsten carbide and carbon with chromium interlayer between substrate and Diamond like Carbon coating.

Plug: Titanium alloy composed of Titanium-6 Aluminium-4 Vanadium ELI (Extra Low Interstitial) Alloy. Detailed chemical composition of the Titanium alloy is Titanium balanced with 5.50 - 6.50 wt.% Aluminium, 3.5-4.5 wt.% Vanadium, max. 0.25 wt.% Iron, max. 0.13 wt.% Oxygen, max. 0.08 wt.% Carbon, max. 0.05 wt.% Nitrogen, max. 0.012 wt.% Hydrogen, (max. -maximum value).

Cautions

General

Treatment by means of implants may lead to loss of bone, biologic, or mechanical failures including fatigue fracture of implants.

Close cooperation between surgeon, restorative dentist, and dental laboratory technician is essential for a successful implant treatment.

Cover screws must only be used with compatible Nobel Biocare instruments and components. Use of instruments or components that are not intended to be used in combination with cover screws can lead to product failure, damage to tissue, or unsatisfactory esthetic results.

When using a new device/treatment method for the first time, working with a colleague who is experienced with the new device/treatment method may help avoid possible complications. Nobel Biocare has a global network of mentors available for this purpose.

Before Surgery

Careful psychological and physiological evaluation followed by clinical and radiological examination must be performed on the patient prior to surgery to determine the suitability of the patient for treatment.

Special attention must be given to patients who have local or systemic factors that could interfere with the healing process of either bone or soft tissue or the osseointegration process (e.g. cigarette smoking, poor oral hygiene, uncontrolled diabetes, oro-facial radiotherapy, steroid therapy, and infections in the neighboring bone). Special caution is advised in patients who receive bisphosphonate therapy.

In general, implant placement and prosthetic design must accommodate individual patient conditions. In case of bruxism, other parafunctional habits or unfavorable jaw relationships, reappraisal of the treatment option may be considered.

The device has not been evaluated in pediatric/adolescent patients and is not recommended for use in children. Routine treatment is not recommended until the end of the juvenile jaw bone growth phase has been properly documented.

Pre-operative hard tissue or soft tissue deficits may yield a compromised esthetic result or unfavorable implant angulations.

All components, instruments, and tooling used during the clinical and/or laboratory procedure must be maintained in good condition and care must be taken that instrumentation does not damage implants or other components.

At Surgery

Care and maintenance of sterile instruments are crucial for successful treatment. Sterilized instruments not only safeguard patients and staff against infection but are also essential for the outcome of the total treatment.

Because of the small sizes of the devices, care must be taken that they are not swallowed or aspirated by the patient. It is appropriate to use specific supporting tools to prevent aspiration of loose parts (e.g. gauze, a dental dam, or throat shields).

After the implant placement, the surgeon's evaluation of bone quality and primary stability will determine when the implants may be loaded. Lack of adequate quantity and/or quality of remaining bone, infection, and generalized diseases may be potential causes for failure of osseointegration both immediately after surgery, or after osseointegration is initially achieved.

After Surgery

To help ensure a successful long term-treatment outcome, it is advised to provide comprehensive regular patient follow up after implant treatment and to inform the patient about appropriate oral hygiene.

Intended Users and Patient Groups

Cover Screws are to be used by dental health care professionals.

Cover Screws are to be used in patients subject to dental implant treatment.

Clinical Benefits and Undesirable Side Effects

Clinical Benefits Associated with Devices in the IFU

Cover Screws are a component of treatment with a dental implant system and/or dental crowns and bridges. As a clinical benefit of treatment, patients can expect to have their missing teeth replaced and/or crowns restored.

Undesirable Side Effects Associated with Cover Screws

During placement of this device the pharyngeal (gag) reflex may be triggered in patients with a sensitive gag reflex. During the submerged healing period, bone may grow over the cover screw. In some cases, cover screws may get exposed prematurely.

Where required per the European Medical Device Regulation (MDR; EU 2017/745) a Summary of Safety and Clinical Performance document (SSCP) is available for the cover screws.

The SSCP can be obtained at the following website:
ec.europa.eu/tools/eudamed¹

¹ Website available upon launch of the European Database on Medical Devices (EUDAMED)

Notice regarding serious incidents

For a patient/user/third party in the European Union and in countries with an identical regulatory regime (Regulation 2017/745/EU on Medical Devices); if, during the use of this device or as a result of its use, a serious incident has occurred, please report it to the manufacturer and to your national authority. The contact information for the manufacturer of this device to report a serious incident is as follows:

Nobel Biocare AB
www.nobelbiocare.com/complaint-form

Surgical Procedure and Handling Procedure

1. Select the appropriate Cover Screw based on the implant connection and platform type.
2. Connect the Cover Screw to the implant and hand tighten using a Unigrip™ Screwdriver, a Driver Brånemark System® Hexagon or an Omnigrip™ Mini, depending on the connection type (see Figure 1).

Note For Brånemark System® cover screw a hex screwdriver should be used for tightening.

Caution Tighten the Cover Screw only finger-tight to avoid excessive loads that might damage the Cover Screw parts.

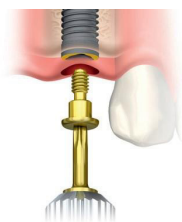


Figure A – Tightening the Cover Screw

3. To remove the Cover Screw, untighten the screw by hand using the appropriate screwdriver.
4. In the event the Cover Screw becomes stuck or is broken, abutment screw retrieval instruments can be used. Refer to Nobel Biocare publication Compatibility Information by navigating to ifu.nobelbiocare.com for the compatible instruments and Nobel Biocare IFU1043 for further details.

Table 1 – Compatible Abutment Screw Retrieval Instruments

| Cover Screw | Abutment Screw Retrieval Instruments |
|--------------------------------------|---|
| Cover Screw Brånemark System® NP | Abutment Screw Retrieval Reverse Drill 3.0/NP Abutment Screw Retrieval Instrument 3.0/NP |
| Cover Screw Brånemark System® RP | Abutment Screw Retrieval Reverse Drill RP/WP/6.0 Abutment Screw Retrieval Instrument RP/WP/6.0 |
| Cover Screw Brånemark System® WP | Abutment Screw Retrieval Reverse Drill RP/WP/6.0 Abutment Screw Retrieval Instrument RP/WP/6.0 |
| Cover Screw Conical Connection 3.0 | Abutment Screw Retrieval Reverse Drill CC 3.0/NP & TCC NP/RP Abutment Screw Retrieval Instrument CC 3.0/NP & TCC NP/RP |
| Cover Screw Conical Connection NP | Abutment Screw Retrieval Reverse Drill CC 3.0/NP & TCC NP/RP Abutment Screw Retrieval Instrument CC 3.0/NP & TCC NP/RP |
| Cover Screw Conical Connection RP | Abutment Screw Retrieval Reverse Drill RP/WP/6.0 Abutment Screw Retrieval Instrument RP/WP/6.0 |
| Cover Screw Conical Connection WP | Abutment Screw Retrieval Reverse Drill RP/WP/6.0 Abutment Screw Retrieval Instrument RP/WP/6.0 |
| Cover Screw NobelReplace® NP | Abutment Screw Retrieval Reverse Drill 3.0/NP Abutment Screw Retrieval Instrument 3.0/NP |
| Cover Screw NobelReplace® RP | Abutment Screw Retrieval Reverse Drill RP/WP/6.0 Abutment Screw Retrieval Instrument RP/WP/6.0 |
| Cover Screw NobelReplace® WP | Abutment Screw Retrieval Reverse Drill RP/WP/6.0 Abutment Screw Retrieval Instrument RP/WP/6.0 |
| Cover Screw NobelReplace 6.0 | Abutment Screw Retrieval Reverse Drill RP/WP/6.0 Abutment Screw Retrieval Instrument RP/WP/6.0 |
| Cover Screw Nobel Biocare N1™ TCC NP | Abutment Screw Retrieval Reverse Drill CC 3.0/NP & TCC NP/RP Abutment Screw Retrieval Instrument CC 3.0/NP & TCC NP/RP |
| Cover Screw Nobel Biocare N1™ TCC RP | Abutment Screw Retrieval Reverse Drill CC 3.0/NP & TCC NP/RP Abutment Screw Retrieval Instrument CC 3.0/NP & TCC NP/RP |

Sterility and Reusability Information

Cover screws have been sterilized using irradiation and are intended for single use only. Do not use after the labeled expiration date.

Warning Do not use device if the packaging has been damaged or previously opened as the device sterility and/or integrity may be compromised.

Caution Cover Screws are single use products and must not be reprocessed. Reprocessing could cause loss of mechanical, chemical and/or biological characteristics. Reuse could cause local or systemic infection.

Magnetic Resonance (MR) Safety Information

These products are fabricated from a metal material which can be affected by MR energy. For further information refer to Nobel Biocare publication **Magnetic Resonance (MR) Safety Information** by navigating to ifu.nobelbiocare.com.

Performance Requirements and Limitations

To achieve the desired performance, the devices must only be used with the products described in this Instructions for Use and/or in the Instructions for Use for other compatible Nobel Biocare products, and in accordance with the Intended Use for each product. To confirm the compatibility of products which are intended to be used in conjunction with the devices, check the color coding, dimensions, lengths, connection type and/or any direct marking as applicable on the products or product labeling.

Facilities and Training

It is strongly recommended that new and experienced users of Nobel Biocare products always go through special training before using a new product for the first time. Nobel Biocare offers a wide range of courses for various levels of knowledge and experience. For more information please visit www.nobelbiocare.com.

Storage, Handling and Transportation

The device must be stored and transported in dry conditions in the original packaging at room temperature and not exposed to direct sunlight. Incorrect storage and transportation may influence device characteristics leading to failure.

Disposal

Safely discard potentially contaminated or no longer usable medical devices as healthcare (clinical) waste in accordance with local healthcare guidelines, country and government legislation or policy.

Separation, re-cycling or disposal of packaging material shall follow local country and government legislation on packaging and packaging waste, where applicable.

Manufacturer and Distributor Information

| | |
|--|--|
| Manufacturer  | Nobel Biocare AB PO Box 5190, 402 26 Västra Hamngatan 1 Göteborg 411 17 Sweden www.nobelbiocare.com |
| UK Responsible Person <div>UK RP</div> | Nobel Biocare UK Ltd 4 Longwalk Road Stockley Park Uxbridge UB11 1FE United Kingdom |
| Distributed in Turkey by | EOT Dental Sağlık Ürünleri ve Dış Ticaret A.Ş. Nispetiye Mah. Aytar Cad. Metro İş Merkezi No: 10/7 Beşiktaş İSTANBUL Phone: +90 2123614901, Fax: +90 2123614904 |
| CE Mark for Class IIb Devices |  |
| UKCA Mark for Class IIb Devices |  |

Note Refer to the product label to determine the applicable conformity marking for each device.

Note Regarding Canadian device licensure, not all products described in the IFU may have a device licence according to Canadian law.

Basic UDI-DI Information

| Product | Basic UDI-DI Number |
|--------------------------------------|---------------------|
| Cover Screw Brånemark System® NP | 73327470000001336W |
| Cover Screw Brånemark System® RP | |
| Cover Screw Brånemark System® WP | |
| Cover Screw Conical Connection 3.0 | 73327470000001316S |
| Cover Screw Conical Connection NP | |
| Cover Screw Conical Connection RP | |
| Cover Screw Conical Connection WP | |
| Cover Screw Nobel Biocare N1™ TCC NP | 73327470000002116R |
| Cover Screw Nobel Biocare N1™ TCC RP | |
| Cover Screw NobelReplace® | 73327470000002226W |

Legal Statements

EN All rights reserved.

Nobel Biocare, the Nobel Biocare logotype and all other trademarks used in this document are, if nothing else is stated or is evident from the context in a certain case, trademarks of Nobel Biocare. Product images in this folder are not necessarily to scale. All product images are for illustration purposes only and may not be an exact representation of the product.

Symbols Glossary

Please refer to the packaging label for the applicable symbols related to the product. On the packaging label you may encounter various symbols to convey a specific information about the product and/or its use. For further information refer to Nobel Biocare publication to the **Symbols Glossary** by navigating to ifu.nobelbiocare.com.