

Ball Abutment

Instructions for use



Important: Please read.

Disclaimer of liability:

This product is part of an overall concept and may only be used in conjunction with the associated original products according to the instructions and recommendation of Nobel Biocare. Non-recommended use of products made by third parties in conjunction with Nobel Biocare products will void any warranty or other obligation, express or implied, of Nobel Biocare. The user of Nobel Biocare products has the duty to determine whether or not any product is suitable for the particular patient and circumstances. Nobel Biocare disclaims any liability, express or implied, and shall have no responsibility for any direct, indirect, punitive or other damages, arising out of or in connection with any errors in professional judgment or practice in the use of Nobel Biocare products. The user is also obliged to study the latest developments in regard to this Nobel Biocare product and its applications regularly. In cases of doubt, the user has to contact Nobel Biocare. Since the utilization of this product is under the control of the user, it is his/her responsibility. Nobel Biocare does not assume any liability whatsoever for damage arising thereof. Please note that some products detailed in this Instruction for Use may not be regulatory cleared, released or licensed for sale in all markets.

Description:

A premanufactured dental implant abutment directly connected to the endosseous dental implant intended for use as an aid in prosthetic rehabilitation. An abutment together with an adjustable Goldcap to be used for full-arch restorations.

Internal tri-channel connection for: NobelReplace®, Replace Select™ and NobelSpeedy® Replace.

External hex connection for: Brånemark System® and NobelSpeedy® Groovy.

Intended use:

Dental implant abutments are intended to be used in the upper or lower jaw and used for supporting tooth replacements to restore chewing function.

The abutments in combination with endosseous implants are used as the foundation for anchoring tooth replacements in either jaw.

Indications:

For fully edentulous jaws retaining a tissue-supported overdenture.

Indicated for patient with extensive bone or soft tissue loss, with compromised manual dexterity or phonetic concerns. Also to be used when implants are placed too posterior to allow retention by means of a bar. Allows for misangulation up to 30° between implants.

Contraindications:

Ball Abutment is contraindicated for patients:

- who are medically unfit for an oral surgical procedure.
- in whom adequate sizes, numbers or desirable positions of implants are not reachable to achieve safe support of functional or eventually parafunctional loads.
- allergic or hypersensitiv to titanium alloy Ti-6Al-4V (titanium, aluminum, vanadium) or gold alloy (gold, platinum, palladium, iridium).

Cautions:

Pre-operative hard tissue or soft tissue deficits may yield a compromised esthetic result or unfavorable implant angulations.

To secure the long term treatment outcome it is advised to provide comprehensive regular patient follow up after implant treatment and to inform about appropriate oral hygiene.

All instruments and tooling used in surgery must be maintained in good condition and care must be taken that instrumentation does not damage implants or other components.

Special attention has to be given to patients who have local or systemic factors that could interfere with the healing process of either bone or soft tissue or the osseointegration process (e.g., cigarette smoking, poor oral hygiene, uncontrolled diabetes, oro-facial radiotherapy, steroid therapy, infections in the neighboring bone).

Special caution is advised in patients who receive bisphosphonate therapy.

In general, implant placement and prosthetic design must accommodate individual patient conditions. In case of bruxism or unfavorable jaw relationships reappraisal of the treatment option may be considered.

Never exceed **15Ncm** prosthetic tightening torque for the abutment. Overtightening of abutment may lead to a screw fracture.

Because of the small size of the devices, care must be taken that they are not swallowed or aspirated by the patient.

It is strongly recommended that clinicians, new as well as experienced implant users, always go through special training before undertaking a new treatment method.

Nobel Biocare offers a wide range of courses for various levels of knowledge and experience. For more info please visit www.nobelbiocare.com.


Working the first time with a colleague, experienced with the new device/treatment method, avoids eventual complications. Nobel Biocare has a global network of mentors available for this purpose.

Close cooperation between surgeon, restorative dentist and dental laboratory technician is essential for a successful implant treatment.


Handling procedure:

Clinical procedure

1. Choose appropriate abutment height and connect the abutments to implants (A).

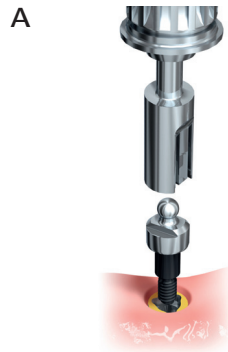


2. Tighten abutments to **15Ncm** using the Manual Torque Wrench Prosthetic and Screwdriver Machine Ball Abutment (B). If the Adapter for NobelReplace® WP to RP is used, push the ball abutment RP through the adapter into the implant and tighten to **15 Ncm** as described above (C).

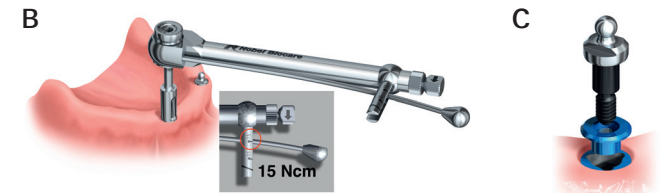

3. Take an impression using a custom tray (D).
4. Position the ball abutment replicas into the impression.

Laboratory procedure

5. Fabricate a model and produce an occlusal rim (E).

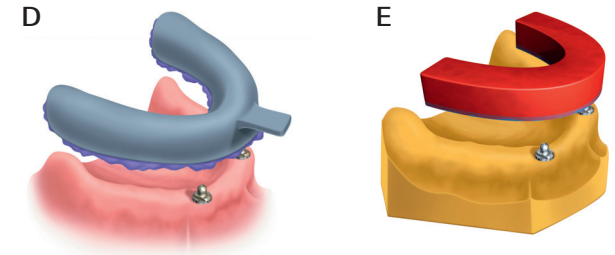

Clinical procedure


6. Check jaw relationship by using the occlusal rim. Make sure the rim does not make contact with ball abutments. Record the jaw relationship.

Laboratory procedure

7. Produce a tooth set-up in wax (F).

Clinical procedure


8. Evaluate wax set-up (G). Make sure that the rim does not make contact with the ball abutments.

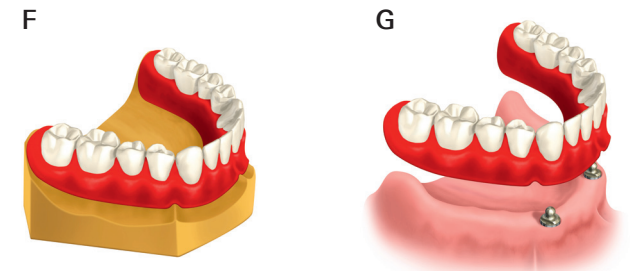

Laboratory procedure

9. Finalize the overdenture and process the gold caps into the overdenture (H).

H

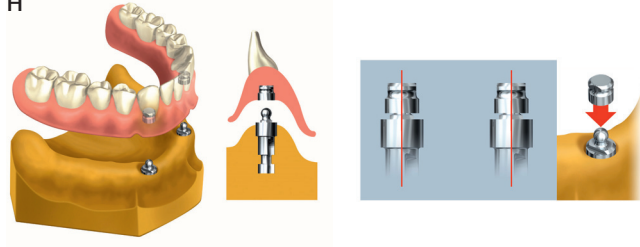


Clinical procedure

10. Use the Screwdriver Machined Ball Abutment and Manual Torque Wrench Prosthetic to verify the tightening of the abutment to **15Ncm**.

Caution: Never exceed recommended maximum **15Ncm** prosthetic tightening torque for the abutment screw. Over tightening of abutment may lead to a screw fracture.

11. Place the overdenture and verify occlusion. If applicable, adjust the retention and confirm the hinge axis movement. Retention force of Gold Caps should be adjusted to desired retention by turning the lamellae retention insert clockwise (increasing) or counter-clockwise (decreasing) using the Screwdriver/Activator for Gold Cap (I).

I

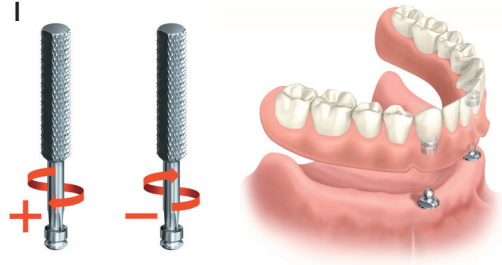

Note: Do not turn the Screwdriver/Activator more than one turn.

Materials:

Ball Abutment: Titanium alloy 90% Ti 6% Al 4% V.

Gold Cap: Gold alloy 60% Au, 19% Pt, 20% Pd, 1% Ir.

Cleaning and sterilization instructions:

Ball Abutment is delivered sterile for single use only prior to the labeled expiration date.

Warning: Do not use if package is damaged or previously opened.

Caution: Ball Abutment is a single use product that must not be reprocessed. Reprocessing could cause loss of mechanical, chemical and / or biological characteristics. Reuse could cause cross contamination.

MRI safety information:

Please note that the product has not been evaluated for safety and compatibility in the MR environment. The product has not been tested for heating or migration in the MR environment.

For additional information on Magnetic Resonance Imaging, please consult the "Cleaning & Sterilization Guidelines including MRI Information of Nobel Biocare Products" available at www.nobelbiocare.com or request latest printed version from a Nobel Biocare representative.

Storage and handling:

The product must be stored in a dry place in the original packaging at room temperature and not exposed to direct sunlight. Incorrect storage may influence device characteristics leading to failure.

Disposal:

Disposal of the device shall follow local regulations and environmental requirements, taking different contamination levels into account.



Manufacturer: Nobel Biocare AB, Box 5190, 402 26

Västra Hamngatan 1, 411 17 Göteborg, Sweden.

Phone: +46 31 81 88 00. Fax: +46 31 16 31 52. www.nobelbiocare.com

Canada license exemption: Please note that not all products may have been licensed in accordance with Canadian law.

Prescription device: Rx only

Caution: Federal law restricts this device to sale by or on the order of a licensed physician or dentist.

CE 0086

Rx Only

STERILE R

Sterile using irradiation



Consult instructions for use



Use-by date



Do not re-use

LOT

Batch code



Do not use if package is damaged

EN All rights reserved.

Nobel Biocare, the Nobel Biocare logotype and all other trademarks used in this document are, if nothing else is stated or is evident from the context in a certain case, trademarks of Nobel Biocare. Product images in this folder are not necessarily to scale.

IFU1024 000 00