

Guided (Pilot Drill) Sleeves, Guided Anchor Pin Sleeves, Mounting Tools



Important – Disclaimer of Liability

This product is part of an overall concept and may only be used in conjunction with the associated original products according to the instructions and recommendation of Nobel Biocare. Non-recommended use of products made by third parties in conjunction with Nobel Biocare products will void any warranty or other obligation, express or implied, of Nobel Biocare. The user of Nobel Biocare products has the duty to determine whether or not any product is suitable for the particular patient and circumstances. Nobel Biocare disclaims any liability, express or implied, and shall have no responsibility for any direct, indirect, punitive or other damages, arising out of or in connection with any errors in professional judgment or practice in the use of Nobel Biocare products. The user is also obliged to study the latest developments in regard to this Nobel Biocare product and its applications regularly. In cases of doubt, the user has to contact Nobel Biocare. Since the utilization of this product is under the control of the user, they are his/her responsibility. Nobel Biocare does not assume any liability whatsoever for damage arising thereof.

Please note that some products detailed in this Instruction for Use may not be regulatory cleared, released or licensed for sale in all markets.

Description

Guided (Pilot Drill) Sleeves and Guided Anchor Pin Sleeves

The Guided (Pilot Drill) Sleeves are highly precise cylinders embedded in dental surgical templates to define the position, direction and height/depth of surgical sites. The surgical template is a plastic device that is to be fitted onto the soft tissue of the patient. In partial and single cases, it also fits onto the remaining teeth of the jaw. The Guided (Pilot Drill) Sleeves have the long axis identical to the planned long axis of the implant. The level of the outer shoulder of the Guided (Pilot Drill) Sleeve defines the depth of the surgical preparation and implant position since there is a pre-defined relationship between this level and the implant/abutment interface.

Also, embedded in the surgical template are Guided Anchor Pin Sleeves. These sleeves guide the preparation and installation of anchor pins. The anchor pins are thin rods of metal positioned close to horizontally into the jawbone to secure the surgical template in its intended position during implant surgery. Guided Anchor Pin Sleeves are optional.

Mounting Tools for Guided (Pilot Drill) Sleeves

In decentralized/local production of surgical templates, the Mounting Tools for Guided (Pilot Drill) Sleeves can be used to glue the Guided (Pilot Drill) Sleeves into the printed/milled/stereolithography (SLA) produced surgical template by the laboratory or dentist. The Mounting Tools for Guided (Pilot Drill) Sleeves will position the Guided (Pilot Drill) Sleeves accurately flush with the top shoulder of the surgical template's sleeve support material. For fixation of the sleeves a biocompatible glue/cement/adhesive agent is needed.

The Mounting Tool Pins for Guided Pilot Drill Sleeves have either a 1.5 mm or 2.0 mm diameter and the Mounting Tool Base is designed to be used with either diameter.

Due to their design, the Mounting Tools for Guided Pilot Drill Sleeves cannot be used in combination of pre-fabrication of master casts.

The Mounting Tools for Guided Sleeves (guided cylinder with pin, implant replica) can be used in a dental lab procedure when making a stone model based on the surgical template. The guided cylinders have an outer diameter that matches the internal dimensions of the Guided Sleeves (NP, RP, and 6.0/WP) and are inserted into the Guided Sleeves, embedded in the surgical template, and an implant replica from the Nobel Biocare standard assortment is screwed onto the pin that goes through the guided cylinder. Since there is a well-defined seating of the guided cylinder in the Guided Sleeve, the head of the implant replica is positioned in the same position in relationship to the surgical template, as the implant will be installed into the bone of the patient. This enables a stone model to be cast using the underside of the surgical template, corresponding to the soft tissue surface of the patient, and the attached implant replicas, corresponding to the pre-defined sites of the implants. The implant connection portion of the guided cylinder should match the available implant connections of tri-channel, external hex and conical connection.

In the analog workflow, the product is attached together with a Guided Sleeve to the implants installed in the die cast model. In doing so, the Guided Sleeves are positioned in a correct position and ready to be molded into a lab-made surgical template.

Nobel Biocare products are intended and available to be used in a variety of configurations. For further information refer to Nobel Biocare publication Compatibility Information by navigating to ifu.nobelbiocare.com.

Intended Use/Intended Purpose

- The Guided (Pilot Drill) Sleeves are intended to be used in combination with surgical templates to facilitate implant installation by guiding the required instruments in the planned direction.
- The Guided Pilot Drill Sleeves are intended to guide the first drill (either 1.5 mm or 2.0 mm).
- The Mounting Tools are intended to be used in the dental laboratory procedure to enable the correct positioning of the Guided (Pilot Drill) Sleeves in the surgical template.
- The Guided Anchor Pin Sleeves are intended to be used in surgical templates to establish secure fixation and stability of the surgical template, by guiding the preparation and installation of anchor pins.

Indications

The components (Guided (Pilot Drill) Sleeves, Guided Anchor Pin Sleeves and Mounting Tools) are indicated for the assembly within a surgical template for the treatment of edentulous and partially edentulous jaws (including patients missing a single tooth) for placement of dental implants, if indicated in combination with immediate function to restore esthetics and functionality (e.g. masticatory, speech). The following prerequisites must be fulfilled:

- Adequate amount of jawbone volume (height and width).
- Quality of jawbone must be judged as adequate.
- Adequate mouth opening (minimum 40 mm) to accommodate guided surgery tooling.

- Exclusion of compromised diseases in conflict with dental implant treatment.
- Adequate compliance.

Contraindications

It is contraindicated to use a surgical template in:

- Patients who are medically unfit for an oral surgical procedure.
- Patients in whom adequate sizes, numbers or desirable positions of implants are not reachable to achieve safe support of functional or eventually parafunctional loads.
- Patients who are allergic or hypersensitivity to commercially pure titanium grade 4, stainless steel, surgical template material or bonding material.

There are no specific contraindications related to the use of the Mounting Tools.

Materials

The Guided (Pilot Drill) Sleeves and Guided Anchor Pin Sleeves are made from stainless steel. The Mounting Tools for the Guided Pilot Drill Sleeves are made of stainless steel. The Mounting Tools for the Guided Sleeves are made of titanium.

Recommendation for materials

For surgical template fabrication

The user must use a certified material, intended for printing surgical guides and follow the manufacturer's instructions for use, as well as the recommended parameters and procedures.

The minimum mechanical properties for the surgical template material are provided in Table 1. In addition, the material must be biocompatible.

Table 1 – Minimum mechanical properties of surgical template material

Post-cured Data	Metric	Method
Tensile Strength	≥41 MPa	ASTM D 638
Tensile Modulus	≥2030 MPa	ASTM D 638
Elongation at Break	4 – 7 %	ASTM D 638
Heat Deflection Temperature		ASTM D 638
at 66Psi	46°C	
at 264Psi	41°C	
Flexural Strength	≥50 MPa	ISO 20795-1/ASTM D 790
Flexural Modulus	≥1500 MPa	ISO 20795-1/ASTM D 790
Hardness Shore	≥80 D	ASTM D2240

For sleeves bonding into the surgical template

The minimum requirements for the mechanical properties of the glue to be used for bonding the sleeves into the surgical template are provided in Table 2.

Table 2 – Minimum mechanical properties of bonding material

Mechanical Property	Acceptable Range/Level
Compressive Strength	≥ 200 MPa
Flexural Strength	> 2000 MPa
Water Sorption	≤ 200 µg/mm ³

Additional requirements for the bonding material

- The bonding material should be biocompatible and suitable for dental applications.
- The bonding material should be able to bond metals to polymers.

Cautions

General

One hundred percent implant success cannot be guaranteed. In particular, non-observance of the product's indications for use and the surgical/handling procedure(s) may result in failure.

Treatment by means of implants may lead to loss of bone, biologic or mechanical failures including fatigue fracture of implants.

Close cooperation between surgeon, restorative dentist and dental laboratory technician is essential for a successful implant treatment.

Nobel Biocare Guided Pilot Drill Sleeves, Guided Sleeves, and Guided Anchor Pin Sleeves must only be used with compatible Nobel Biocare instruments and/or components and/or prosthetic components. Use of instruments and/or components and/or prosthetic components that are not intended to be used in combination with Nobel Biocare Guided Pilot Drill Sleeves, Guided Sleeves, and Guided Anchor Pin Sleeves can lead to product failure, damage to tissue, or unsatisfactory esthetic results.

When using a new device/treatment method for the first time, working with a colleague who is experienced with the new device/treatment method may help avoid possible complications. Nobel Biocare has a global network of mentors available for this purpose.

Before Surgery

Careful psychological and physiological evaluation, followed by clinical and radiological examination must be performed on the patient prior to surgery to determine the suitability of the patient for treatment. Special attention must be given to patients who have local or systemic factors that could interfere with the healing process of either bone or soft tissue or the osseointegration process (e.g. cigarette smoking, poor oral hygiene, uncontrolled diabetes, oro-facial radiotherapy, steroid therapy, infections in the neighboring bone). Special caution is advised in patients who receive bisphosphonate therapy.

In general, implant placement and prosthetic design must accommodate individual patient conditions. In case of bruxism, other parafunctional habits or unfavorable jaw relationships, reappraisal of the treatment option may be considered.

Pre-operative hard tissue or soft tissue deficits may yield a compromised esthetic result or unfavorable implant angulations.

All components, instruments and tooling used during the clinical and/or laboratory procedure must be maintained in good condition and care must be taken that instrumentation does not damage implants or other components.

At Surgery

Small diameter implants and angled abutments are not recommended for the posterior region.

Care and maintenance of sterile instruments are crucial for a successful treatment. Sterilized instruments not only safeguard your patients and staff against infection but are also essential for the outcome of the total treatment.

Because of the small sizes of the devices, care must be taken that they are not swallowed or aspirated by the patient. It is appropriate to use specific supporting tools to prevent aspiration of loose parts (e.g. gauze, dental dam, or throat shield).

After the implant placement, the surgeon's evaluation of bone quality and primary stability will determine when implants may be loaded. Lack of adequate quantity and/or quality of remaining bone, infection and generalized diseases may be potential causes for failure of osseointegration both immediately after surgery, or after osseointegration is initially achieved.

After Surgery

To help ensure a successful long term-treatment outcome, it is advised to provide comprehensive regular patient follow up after implant treatment and to inform the patient about appropriate oral hygiene.

Intended Users and Patient Groups

- Guided Pilot Drill Sleeves, Guided Sleeves, and the Guided Anchor Pin Sleeve are to be used by laboratory professionals.
- Guided Pilot Drill Sleeves, Guided Sleeves, and the Guided Anchor Pin Sleeve are to be used in patients subject to dental implant treatment.
- Mounting Tools are to be used by laboratory professionals.

Clinical Benefits and Undesirable Side Effects

Clinical Benefits Associated with Guided Pilot Drill Sleeves, Guided Sleeves, and the Guided Anchor Pin Sleeve

Guided Pilot Drill Sleeves, Guided Sleeves, and the Guided Anchor Pin Sleeve are a component of treatment with a dental implant system and/or dental crowns and bridges. As a clinical benefit of treatment, patients can expect to have their missing teeth replaced and/or crowns restored.

Clinical Benefits Associated with Mounting Tools and Guided Cylinder with Pin

Mounting Tools and Guided Cylinder with Pin are used to produce dental surgical templates which are components of treatment with a dental implant system and/or dental crowns and bridges. As a clinical benefit of treatment, patients can expect to have their missing teeth replaced and/or crowns restored.

Undesirable Side Effects Associated with Guided Pilot Drill Sleeves, Guided Sleeves, and the Guided Anchor Pin Sleeve

During use of these devices the pharyngeal (gag) reflex may be triggered in patients with a sensitive gag reflex.

Undesirable Side Effects Associated with Mounting Tools and Guided Cylinder with Pin

None known.

Notice regarding serious incidents

For a patient/user/third party in the European Union and in countries with an identical regulatory regime (Regulation 2017/745/EU on Medical Devices); if, during the use of this device or as a result of its use, a serious incident has occurred, please report it to the manufacturer and to your national authority. The contact information for the manufacturer of this device to report a serious incident is as follows:

Nobel Biocare AB
www.nobelbiocare.com/complaint-form

Handling Procedure

- 1. Check the decentralized/locally produced surgical template.
 - Check the sleeve seats for material residues and sharp, protruding edges. Remove or smoothen them if you find any.
 - Check that the locally produced surgical template is manufactured from appropriate material: the material should be biocompatible and mechanically fit for purpose. Recommended material properties are listed in table 1.
 - Verify optimal fit on stone model if applicable and/or in patient’s mouth prior to surgery.
 - 2a. Fixing the Guided Pilot Drill Sleeves into the surgical template.
 - Insert the Guided Pilot Drill Sleeve into the sleeve seat of the surgical template.
- Note** Make sure the flat upper part of the sleeve is on the occlusal surface of the surgical template.
- For the correct fixation of the Guided Pilot Drill Sleeve; use the screw part of the Mounting Tool from below (intaglio/ inside of the surgical template) and the thread part from the top (occlusion/outer surface of the surgical template) (see Figure B). Tighten it by hand.

Note The outer diameter of the Guided Pilot Drill Sleeves is the same.

Please confirm that the correct sleeve size is used for each position by comparing the laser marking of the Guided Pilot Drill Sleeves with the planning overview before inserting it into the template (see Figure A).

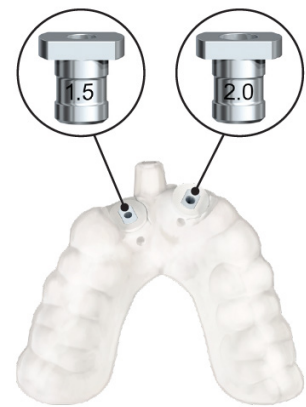


Figure A – Example of a comparison of laser marking on Guided Pilot Drill Sleeves with planning overview

PID	Sleeve name	Article id
1	Guided Pilot Drill Sleeve 1.5 mm	300438
2	Guided Pilot Drill Sleeve 2.0 mm	300440

- 2b. Fixing the Guided Sleeve into the surgical template.
 - Insert the Guided Sleeve into the sleeve seat of the surgical template.
- Note** As the Guided Sleeve is symmetrical, there is no top or bottom.
- For the correct fixation of the Guided Sleeve; use the upper part of the Guided Cylinder with Pin and its lower part to fix the Guided Sleeve. Tighten the lower part of the Guided Cylinder with Pin to the appropriate Implant Replica (see Figure C). Fix all parts together and tighten them by hand or by using a screwdriver (e.g. Unigrip™).
- 2c. Fixing the Guided Anchor Pin Sleeve.
 - Insert the Guided Anchor Pin Sleeve into the sleeve seat of the surgical template.
- Note** As the Guided Anchor Pin Sleeve is symmetrical, there is no top or bottom.
- Make sure that the upper part of the sleeve is set flush into the seat socket of the surgical template (see Figure D).

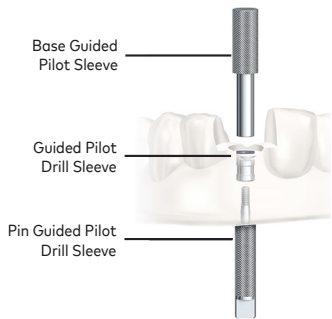


Figure B – Fixation of Guided Pilot Drill Sleeve

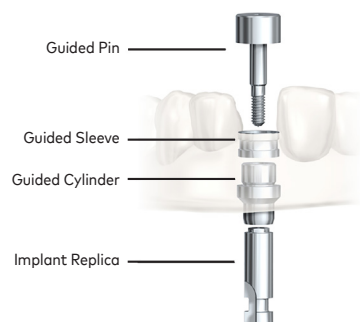


Figure C – Fixation of Guided Sleeve



Figure D – Insertion of Anchor Pin Sleeve

3. Bonding process for the Guided (Pilot Drill) Sleeves and Guided Anchor Pin Sleeves.
 - Visually verify that the sleeves are flush with the top surrounding surface of the surgical template. If they are not flush, remove material as required.

For permanent fixation of the sleeves a biocompatible glue/cement/adhesive agent is needed. The user must use a biocompatible material and follow the manufacturer's instructions for use. Recommended material properties are listed in Table 2.

Note Once all the Guided (Pilot Drill) Sleeve(s) and Guided Anchor Pin Sleeve(s) are in place, glue them into the surgical template.

- Place the tip of the mixing cannula inside the glue hole located in the surgical template.
- Slowly squeeze the bonding material into the glue hole until it completely goes all the way around the Guided (Pilot Drill) Sleeve(s) or Guided Anchor Pin Sleeve(s).
- Repeat the application of bonding material to each Guided (Pilot Drill) Sleeve(s) and Guided Anchor Pin Sleeve(s).
- Once the bonding material has hardened, unscrew the Mounting Tools.
- Visually verify that there is no bonding material present on the top or bottom of the Guided (Pilot Drill) Sleeve(s) or Guided Anchor Pin Sleeve(s).

Dual curing is advised to control the timing of the curing process and ensure material is fully cured.

Caution Introduce only as much bonding material to cover the outer diameter of the Guided (Pilot Drill) Sleeve(s) or Guided Anchor Pin Sleeve(s) in the glue channel. While introducing the bonding material, you should be able to see the progression of it. Avoid to introduce too much bonding material. Leaked bonding material should be removed immediately using a suitable instrument.

Warning Do not introduce bonding material into the internal diameter of the Guided (Pilot Drill) Sleeve(s) or Guided Anchor Pin Sleeve(s). This could have a negative effect on the guided surgery.

For additional information on surgical procedures please consult the NobelGuide Concept Manual available at www.nobelbiocare.com or request the latest printed version from a Nobel Biocare representative.

Sterility and Reusability Information

Delivered non-sterile for single use

Guided (Pilot Drill) Sleeves and Guided Anchor Pin Sleeves are delivered non-sterile for single use only and do not need to be sterilized prior to use.

Caution Guided (Pilot Drill) Sleeves and Guided Anchor Pin Sleeves are single use products and must not be reprocessed. Reprocessing could cause loss of mechanical, chemical and/or biological characteristics. Re-use could cause cross contamination.

The Mounting Tools are used in the dental laboratory only (no intraoral use) and has/have no cleaning and/or sterilization requirements.

Surgical templates must be cleaned and disinfected prior to intraoral use following the procedures in the Cleaning and Disinfection Instructions. During processing in the dental laboratory, templates can be cleaned as necessary without disinfection.

Cleaning and Sterilization Instructions

These products are intended to be cleaned and sterilized. For further information refer to Nobel Biocare publication **Cleaning and Sterilization Instructions** by navigating to ifu.nobelbiocare.com.

Performance Requirements and Limitations

To achieve the desired performance, the devices must only be used with the products described in this Instructions for Use and/or in the Instructions for Use for other compatible Nobel Biocare products, and in accordance with the Intended Use for each product. To confirm the compatibility of products which are intended to be used in conjunction with the devices, check the color coding, dimensions, lengths, connection type and/or any direct marking as applicable on the products or product labeling.

Facilities and Training

It is strongly recommended that new and experienced users of Nobel Biocare products always go through special training before using a new product for the first time. Nobel Biocare offers a wide range of courses for various levels of knowledge and experience. For more information please visit www.nobelbiocare.com.

Storage, Handling and Transportation

The device must be stored and transported in dry conditions in the original packaging at room temperature and not exposed to direct sunlight. Incorrect storage and transportation may influence device characteristics leading to failure.

Disposal

Safely discard potentially contaminated or no longer usable medical devices as healthcare (clinical) waste in accordance with local healthcare guidelines, country and government legislation or policy.

Separation, re-cycling or disposal of packaging material shall follow local country and government legislation on packaging and packaging waste, where applicable.

Manufacturer and Distributor Information

Manufacturer	
	Nobel Biocare AB PO Box 5190, 402 26 Västra Hamngatan 1 Göteborg 411 17 Sweden www.nobelbiocare.com
UK Responsible Person	
UK RP	Nobel Biocare UK Ltd 4 Longwalk Road Stockley Park Uxbridge UB11 1FE United Kingdom
Distributed in Turkey by	
	EOT Dental Sağlık Ürünleri ve Dış Ticaret A.Ş Nispetiye Mah. Aytaç Cad. Metro İş Merkezi No: 10/7 Beşiktaş İSTANBUL Phone: +90 2123614901, Fax: +90 2123614904
CE Mark for Class I Devices	
	
UKCA Mark for Class I Devices	
	

Note Refer to the product label to determine the applicable conformity marking for each device.

Note Regarding Canadian device licensure, not all products described in the IFU may have a device licence according to Canadian law.

Basic UDI-DI Information

Product	Basic UDI-DI Number
Guided Pilot Drill Sleeves 1.5 mm/2.0 mm	733274700000013572
Guided Sleeves NP/RP/6.0/WP	733274700000013572
Guided Anchor Pin Sleeve 1.5 mm	73327470000001957L
Guided Cylinder w Pin Unigrip BmkSyst NP	0733274700000020874
Guided Cylinder w Pin Unigrip BmkSyst RP	0733274700000020874
Guided Cylinder w Pin Unigrip BmkSyst WP	0733274700000020874
Guided Cylinder w Pin Unigrip NobRpl NP	0733274700000020874
Guided Cylinder w Pin Unigrip NobRpl RP	0733274700000020874
Guided Cylinder w Pin Unigrip NobRpl WP	0733274700000020874
Guided Cylinder w Pin Unigrip NobRpl 6.0	0733274700000020874
Guided Cylinder w Pin Conical Connection NP	0733274700000020874
Guided Cylinder w Pin Conical Connection NP 3.5	0733274700000020874
Guided Cylinder w Pin Conical Connection RP 4.3	0733274700000020874
Guided Cylinder w Pin Conical Connection RP 5.0	0733274700000020874
Guided Cylinder with Pin CC WP 5.5	0733274700000020874
Mounting Tool Pin Guided Pilot Sleeve 1.5mm	0733274700000020874
Mounting Tool Pin Guided Pilot Sleeve 2.0mm	0733274700000020874
Mounting Tool Base Guided Pilot Sleeve 1.5/2.0mm	0733274700000020874

Legal Statements

EN All rights reserved.

Nobel Biocare, the Nobel Biocare logotype and all other trademarks used in this document are, if nothing else is stated or is evident from the context in a certain case, trademarks of Nobel Biocare. Product images in this folder are not necessarily to scale. All product images are for illustration purposes only and may not be an exact representation of the product.

Symbols Glossary

Please refer to the packaging label for the applicable symbols related to the product. On the packaging label you may encounter various symbols to convey a specific information about the product and/or its use. For further information refer to Nobel Biocare publication to the **Symbols Glossary** by navigating to ifu.nobelbiocare.com.