

Temporary Snap Abutment, Engaging/Non-Engaging

Instructions for use



Important: Please read.

Disclaimer of liability:

This product is part of an overall concept and may only be used in conjunction with the associated original products according to the instructions and recommendation of Nobel Biocare. Non-recommended use of products made by third parties in conjunction with Nobel Biocare products will void any warranty or other obligation, express or implied, of Nobel Biocare. The user of Nobel Biocare products has the duty to determine whether or not any product is suitable for the particular patient and circumstances. Nobel Biocare disclaims any liability, express or implied, and shall have no responsibility for any direct, indirect, punitive or other damages, arising out of or in connection with any errors in professional judgment or practice in the use of Nobel Biocare products. The user is also obliged to study the latest developments in regard to this Nobel Biocare product and its applications regularly. In cases of doubt, the user has to contact Nobel Biocare. Since the utilization of this product is under the control of the user, they are his/her responsibility. Nobel Biocare does not assume any liability whatsoever for damage arising thereof. Please note that some products detailed in this Instructions for Use may not be regulatory cleared, released or licensed for sale in all markets.

Description:

A premanufactured dental implant abutment directly connected to the endosseous dental implant intended for use as a temporary aid in prosthetic rehabilitation.

Temporary Snap Abutment Engaging, Temporary Snap Abutment Non-Engaging: Internal conical connection for: NobelActive®, NobelReplace® CC and NobelParallel™ CC.

Intended use:

Temporary Snap Abutments are endosseous dental implant abutments that are designed for single use as a temporary prosthesis during the healing process while the permanent prosthesis is fabricated.

Indications:

The Temporary Snap Abutment is intended to be used to fabricate and support provisional restorations that aid in creating an esthetic emergence through the gingiva during the healing period and prior to final restoration. The Temporary Snap Abutment can be used for cement retained or screw-retained provisional restorations. The abutments can be used for single-unit and multi-unit restorations. Use of the Temporary Snap Abutment is not to exceed one hundred and eighty (180) days.

Contraindications:

Temporary Snap Abutments are contraindicated for patients:

- who are medically unfit for an oral surgical procedure.
- in whom adequate sizes, numbers or desirable positions of implants are not reachable to achieve safe support of functional or eventually parafunctional loads.
- allergic or hypersensitive to commercially pure titanium, titanium alloy Ti-6Al-4V (titanium, aluminum, vanadium).

Cautions:

Pre-operative hard tissue or soft tissue deficits may yield a compromised esthetic result or unfavorable implant angulations. To secure the long term treatment outcome it is advised to provide comprehensive regular patient follow up after implant treatment and to inform about appropriate oral hygiene.

All instruments and tooling used in surgery must be maintained in good condition and care must be taken that instrumentation does not damage implants or other components. Special attention has to be given to patients who have local or systemic factors that could interfere with the healing process of either bone or soft tissue or the osseointegration process (e.g., cigarette smoking, poor oral hygiene, uncontrolled diabetes, oro-facial radiotherapy, steroid therapy, infections in the neighboring bone).

Special caution is advised in patients who receive bisphosphonate therapy. In general, implant placement and prosthetic design must accommodate individual patient conditions. In case of bruxism or unfavorable jaw relationships reappraisal of the treatment option may be considered. Never exceed the recommended maximum prosthetic tightening torque for the abutment screw (see clinical procedure). Overtightening of abutment may lead to a screw fracture.

Because of the small sizes of the devices, care must be taken that they are not swallowed or aspirated by the patient. It is recommended to use a rubber dam in order to prevent inhalation of loose parts.

It is strongly recommended that clinicians, new as well as experienced implant users, always go through special training before undertaking a new treatment method. Nobel Biocare offers a wide range of courses for various levels of knowledge and experience. For more information, please visit www.nobelbiocare.com.


Working the first time with a colleague, experienced with the new device/treatment method, avoids eventual complications. Nobel Biocare has a global network of mentors available for this purpose.

Close cooperation between surgeon, restorative dentist and dental laboratory technician is essential for a successful implant treatment.

Temporary Snap Abutments shall be taken out of occlusion and should not be used for full-arch restorations.

Handling instructions:

Modifications of abutments could be performed using copious water irrigation. Extra-oral modification of abutment is recommended. Use a carborundum disk and carbide bur to shorten the abutment. Minimum height is 4 mm (figure 1). An angular correction of the abutment or any further adjustments should not be conducted.

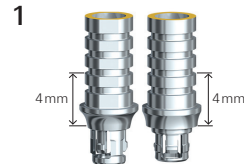

Clinical procedure (Chair-side made provisionals):

1. Connect the Temporary Snap Abutment to the implant and modify the abutment if necessary using copious irrigation.
2. Close the screw access hole.
3. Make a temporary restoration using a pre-fabricated mold with suitable temporary crown and bridge material.
4. Drill a hole through the mold, loosen the screw(s) using a Unigrip™ Screwdriver and remove the restoration.

5. Make final adjustments.
6. Connect the temporary restoration using a Unigrip™ Screwdriver.
7. Tighten abutment to **35 Ncm** using Unigrip™ Machine Screwdriver and Manual Torque Wrench Prosthetic.

Warning: Over tightening of abutment may lead to a screw fracture.

Materials:

Temporary Snap Abutment with Internal conical connection: Titanium alloy 90% Ti, 6% Al, 4%V.

Cleaning and sterilization:

Temporary Snap Abutments are delivered non-sterile for single use and must be cleaned and sterilized prior to use.

Warning: Use of non-sterile components may lead to infection of tissues or infectious diseases.

Warning: Do not use device if the packaging has been damaged or previously opened.

Caution: This is a single use product that must not be reprocessed. Reprocessing could cause loss of mechanical, chemical and/or biological characteristics. Reuse could cause cross contamination.

Manual cleaning, disinfection and drying:

1. Remove debris in lukewarm water and immerse devices in cleaning solution.
2. Scrub with soft bristled nylon brush.
3. Immerse in ultrasonic bath.
4. Rinse and flush with purified or sterile water.
5. Immerse in disinfection solution.
6. Flush internal channels/lumina with disinfection solution.
7. Rinse and flush with purified or sterile water.
8. Dry with compressed air or wipes.

Automated cleaning, disinfection and drying (incl. pre-cleaning):

1. Remove debris in lukewarm water and immerse devices in cleaning solution.
2. Scrub with soft bristled nylon brush.
3. Rinse with tap water.
4. Load devices into thermodisinfector.
5. Run cleaning and disinfection cycle.
6. Run drying cycle.
7. Dry with compressed air or wipes if needed.

Seal single device in a pouch and steam sterilize at 270°F (132°C) for 4 minutes.

Note: FDA cleared sterilization accessories are to be used for the recommended sterilization parameters for wrapping the device supplied non-sterile before user sterilization.

USA

Method	Moist heat sterilization	
	Pre-vacuum	Gravity
Cycle		
Temperature	270°F (132°C)	
Exposure time	4 minutes	15 minutes
Pre-vacuum	3 times < 60 mbar	N/A
Drying time	20-30 minutes	15-30 minutes
Cooling time	10 minutes at room temperature	

Magnetic Resonance (MR) safety information:

The Temporary Snap Abutment has not been evaluated for safety and compatibility in the MR environment. They have not been tested for heating, migration, or image artifact in the MR environment. The safety of the Temporary Snap Abutment in the MR environment is unknown. Scanning a patient who has this device may result in patient injury.

Storage and handling:

The device must be stored and transported in dry conditions in the original packaging at room temperature and not exposed to direct sunlight. Incorrect storage and transportation may influence device characteristics leading to failure.

Disposal:

Disposal of the device shall follow local regulations and environmental requirements, taking different contamination levels into account.



Manufacturer: Nobel Biocare AB, Box 5190, 402 26
Västra Hamngatan 1, 411 17 Göteborg, Sweden.

Phone: +46 31 81 88 00. Fax: +46 31 16 31 52. www.nobelbiocare.com

Canada license exemption: Please note that not all products may have been licensed in accordance with Canadian law.

Prescription device: Rx only

Caution: Federal law restricts this device to sale by or on the order of a licensed physician or dentist.



Rx Only



Consult instructions
for use



Do not re-use



Batch code



Do not use
if package is
damaged

EN All rights reserved.

Nobel Biocare, the Nobel Biocare logotype and all other trademarks used in this document are, if nothing else is stated or is evident from the context in a certain case, trademarks of Nobel Biocare. Product images are not necessarily to scale.